

STERIS Verify Self-Contained Biological Indicators (SCBIs) are the most convenient format of biological indicators. The fully-enclosed indicator is designed to monitor steam and/or ethylene oxide (EO) sterilization cycles. The SCBI system helps eliminate false positives and the need for aseptic handling.



## Verify® Self-Contained Biological Indicators

### Steam and Ethylene Oxide Sterilization

#### FEATURES

Complete system is self-contained, and includes the biological indicator, media and pH indicator

Available in dual species formulation

Dual chemical indicators on label for both steam and EO processes

Reduced incubation time (RIT) studies available

Closed system after the Verify cap is depressed

Vivid color change of pH indicator during incubation to indicate microbial growth

Specifically designed activator set and incubators available

#### BENEFITS

Reduces the chance for false positives due to poor aseptic transfer technique

Reduces inventory

Provides immediate proof of processing

Quicker results, shorter processing times, faster product release and reduced inventory costs

Prevents contamination of vial contents; minimizes media evaporation so a seven day incubation can be used if desired



Allows for easy identification of microbial growth

Easy to process

STERIS SCBIs are designed for validation, requalification or routine use in any size steam or EO sterilizer. Backed by a library of technical documentation, our SCBIs are manufactured in accordance with the International Organization for Standardization (ISO) and the United States Pharmacopeia (USP) standards. Additionally our technical services group is recognized for their technical leadership, expertise and Customer satisfaction. Industry leading products and unparalleled technical service is the STERIS difference.

## VERIFY SELF-CONTAINED BIOLOGICAL INDICATORS

SPECIES	POPULATION	PRODUCT NUMBER	SHELF LIFE	QTY
<i>B. atrophaeus</i> (BA)	6 Log	S3112	12 months	100
<i>G. stearothermophilus</i> (GST)	6 Log	S3111	12 months	100
Dual Species	BA: 6 Log / GST: 5 Log	S3060	12 months	100
Dual Species	BA: 6 Log / GST: 5 Log	S3061	12 months	50

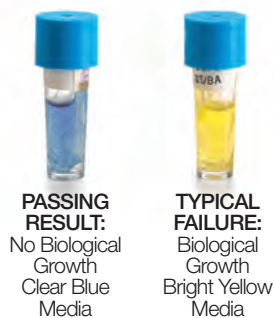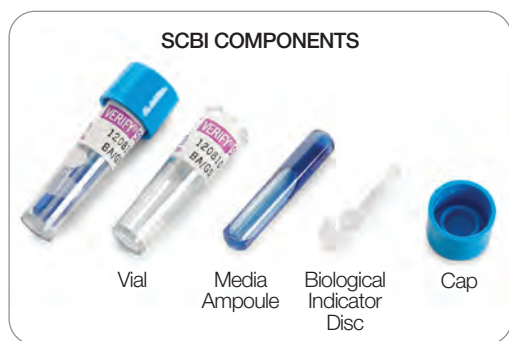

## VERIFY VIAL ACTIVATOR SET

Ensures proper seating of cap and crushing of ampoule during activation step.

DESCRIPTION	PRODUCT NUMBER	QTY
Verify Vial Activator Set	S3075	1



## VERIFY INCUBATORS

Specially designed for use with Verify Self-Contained Biological Indicators. Available in dual temperature model: 35-39°C for EO applications and 55-60°C for steam applications or single temperature model for steam applications.

DESCRIPTION	PRODUCT NUMBER	STANDARD(S)
Dual Incubator (28 well)	S3080	ETL
Steam Incubator (28 well)	S3082	ETL
Steam Incubator (6 well)	S3271	ETL, CSA



STERIS CORPORATION  
5960 Heisley Road  
Mentor, OH 44060, USA  
T +1 800 444 9009

[www.STERISLifeSciences.com](http://www.STERISLifeSciences.com)

### STERIS Offices Worldwide

Asia Pacific	+65 (0) 68 41 7677
Canada	+1 800 444 9009
Finland	+358 9 25851
France	+33 (0) 2 38 84 85 40
Germany	+49 (0) 221 466120
Italy	+39 (0) 2 21303 424
Latin America	+1 800 444 9009
Spain	+34 (0) 916 585 920
United Kingdom	+44 (0) 1256 840400

## STERIS Life Sciences | Science & Solutions for Life

**For over 100 years**, STERIS has been known as a global leader, trusted partner and solutions provider in the field of sterilization and contamination control. **Today**, STERIS continues building on this heritage by advancing the **Science** of sterilization, cleaning and infection control while offering **Solutions** that meet our Customers' needs and high standards. STERIS is dedicated to helping you enhance the **Life** of your patients and the life of your equipment. From Formulated Chemistries, to Capital Equipment, to Parts and Services, STERIS Life Sciences is **Science & Solutions for Life**.

© 2015 STERIS Corporation  
All rights reserved

No part of this publication may be reproduced in any material form without the permission of STERIS Corp.